




nxt™ non-stick bipolar forceps



Contents

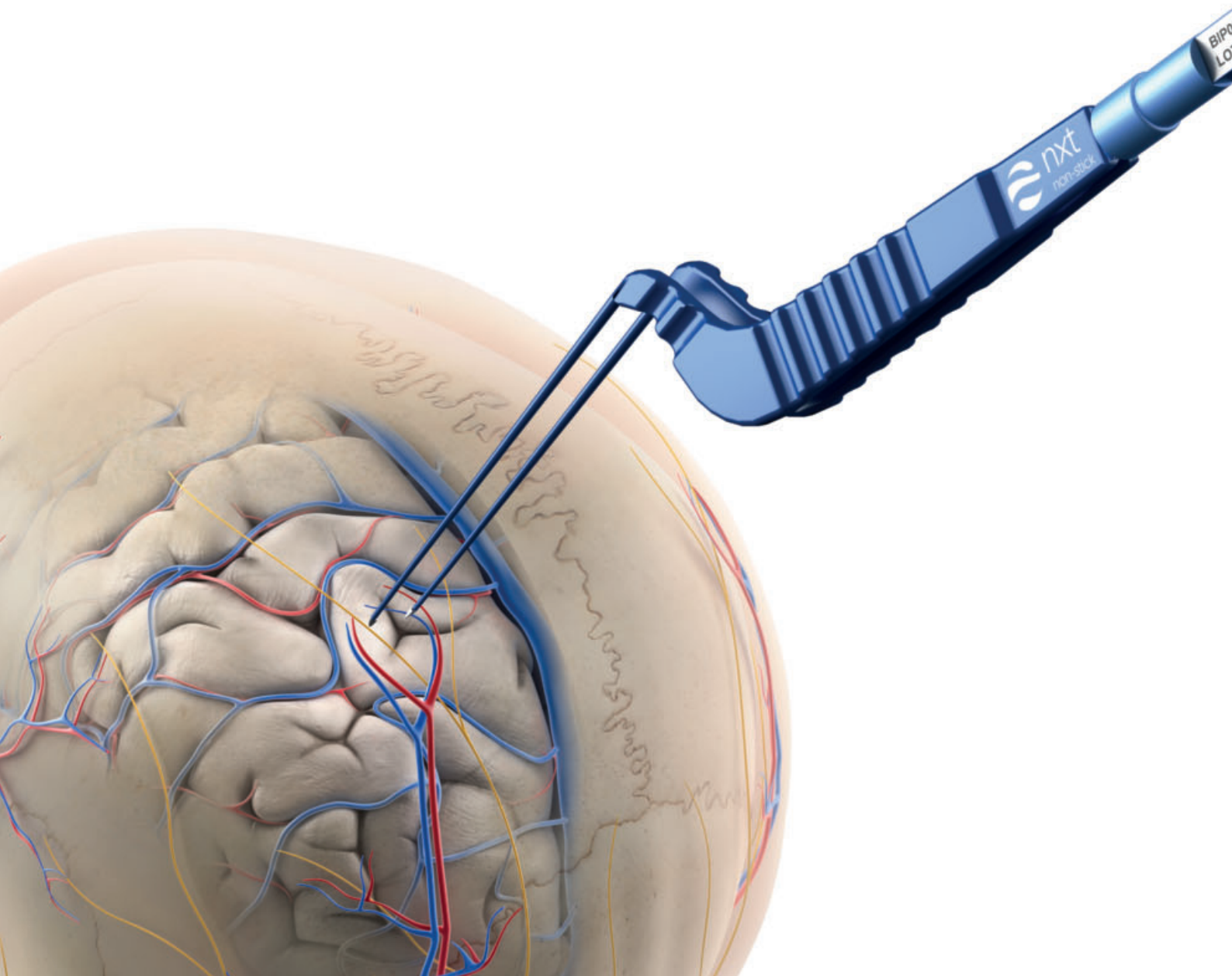

04-28

nxt™ non-stick bipolar forceps
with nxt™ non-stick technology

nxt™ non-stick Bipolarpinzetten
mit nxt™ non-stick Technologie



30-33

nxt™ irrigation non-stick bipolar forceps
with nxt™ irrigation technology

nxt™ non-stick Spülpinzetten
mit nxt™ Spültechnologie



34-39

nxt™ suction bipolar forceps
with nxt™ suction technology

nxt™ Saugpinzetten
mit nxt™ Saugtechnologie

adeor® is a registered trademark of adeor medical AG

All rights of distribution, also by photocopy, reprint or storing and recovery in any kind of data processing device are reserved and require our written approval.

Copyright 2021 adeor medical AG
01-1121-D/E
Printed in Germany
November 2021

Aesculap is a trademark of Aesculap AG
Erbe is a trademark of Erbe Elektromedizin GmbH
Storz is a trademark of Storz Medical AG
Sutter is a trademark of Sutter Medizintechnik GmbH
Martin is a trademark of KLS Martin Group
Valleylab is a trademark of Covidien Surgical Solutions Group
adeor is not affiliated with any of the above.



40-47

Zeppelin™ style round handle bipolar forceps
with nxt™ non-stick technology

Zeppelin™ Rundgriff Bipolarpinzetten
mit nxt™ non-stick Technologie

48-49

bipolar cables

Bipolarkabel



nxt™ non stick - solid silver tips

Superior non-stick technology

For the first time, adeor uses solid silver for its nxt™ non-stick technology. An alloy with unsurpassed thermal and electrical conductivity.

Thanks to these superior thermal conductivity properties, our bipolar instruments offer the industry-leading nxt™ non-stick effect. Another advantage of our forceps is achieved by the polishing process of the silver tips. This further increases the non-stick effect and also facilitates cleaning and reprocessing procedures. The non-stick effect of the adeor bipolar forceps ensures consistent results for the lifetime of the forceps.

nxt™ non stick - Instrumentenspitzen aus reinem Silber


Überlegene Non-Stick Technologie

adeor verwendet für seine nxt™ Non-Stick Technologie erstmals reines Silber, wodurch eine unübertroffene thermische sowie elektrische Leitfähigkeit erreicht wird.

Dieser thermischen Leitfähigkeit von reinem Silber sind die einzigartigen Antihafteigenschaften der nxt™ non-stick Pinzetten zu verdanken. Reines Silber liefert derzeit die besten non-stick Eigenschaften auf dem Markt für Bipolarpinzetten. nxt™ non-stick Pinzetten verfügen im Gegensatz zu Konkurrenzprodukten über polierte Silberspitzen, wodurch der non-stick Effekt noch weiter verstärkt wird. Die glatte Spitzenoberfläche bietet neben dem verbesserten non-stick Effekt weitere Vorteile bei Reinigung und Wiederaufbereitung. Der non-stick Effekt bleibt durch die Verwendung von reinem Silber über die gesamte Lebensdauer der Pinzette erhalten.



nxtTM
 Non-Stick Technology VS. **standard** bipolar forceps

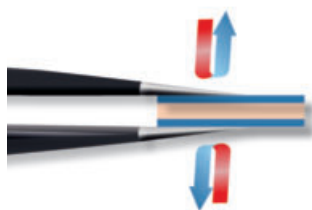

Special nxtTM material

Active heat transport

No sticking of the forceps!

Standard material

Forceps tips overheat

Forceps begin to stick



Rupture of coagulated tissue






nano™ non-stick



nano™ tips (4 mm length)
 nano™ Spitzen (4 mm Länge)

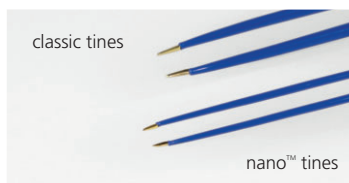

classic tines
 nano™ tines
 classic Ø 3.5 mm vs. nano™ Ø 2 mm tines
 classic Ø 3.5 mm vs. nano™ Ø 2 mm Schenkel



TruPar™ Pin guide



Wide ergonomic grip
 Breiter Ergo-Griff



Completely closed Hygiene+™ connector
 Geschlossener Hygiene+™ Stecker



classic non-stick




TruPar™ Pin guide



Wide ergonomic grip
Breiter Ergo-Griff



Completely closed Hygiene+™ connector
Geschlossener Hygiene+™ Stecker

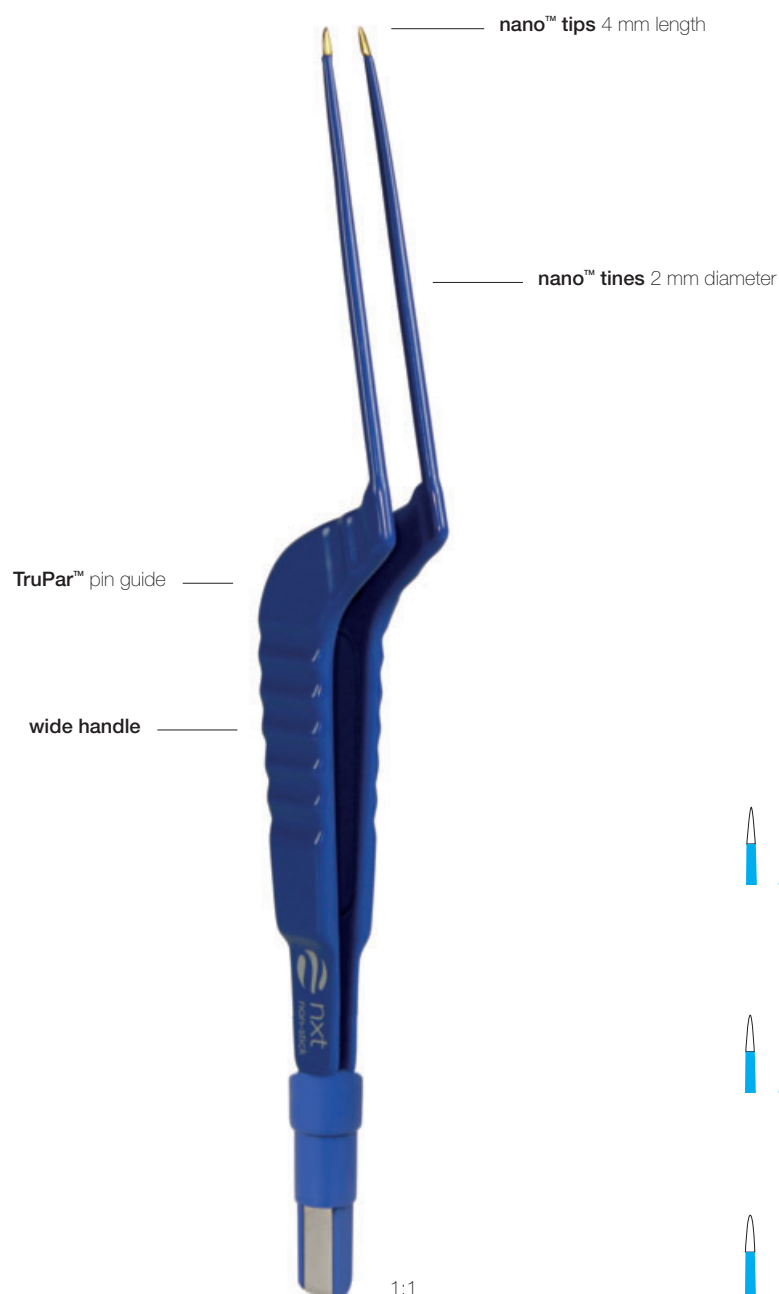

nano™ non-stick



Ø 0.2 x 4 mm nano™ tip straight

nano™ series / working / total length 55 / 170 mm

BIP001EU eu flat connector

BIP001US us pin connector



Ø 0.4 x 4 mm nano™ tip straight

nano™ series / working / total length 55 / 170 mm

BIP002EU eu flat connector

BIP002US us pin connector




Ø 0.7 x 4 mm nano™ tip straight

nano™ series / working / total length 55 / 170 mm

BIP003EU eu flat connector

BIP003US us pin connector

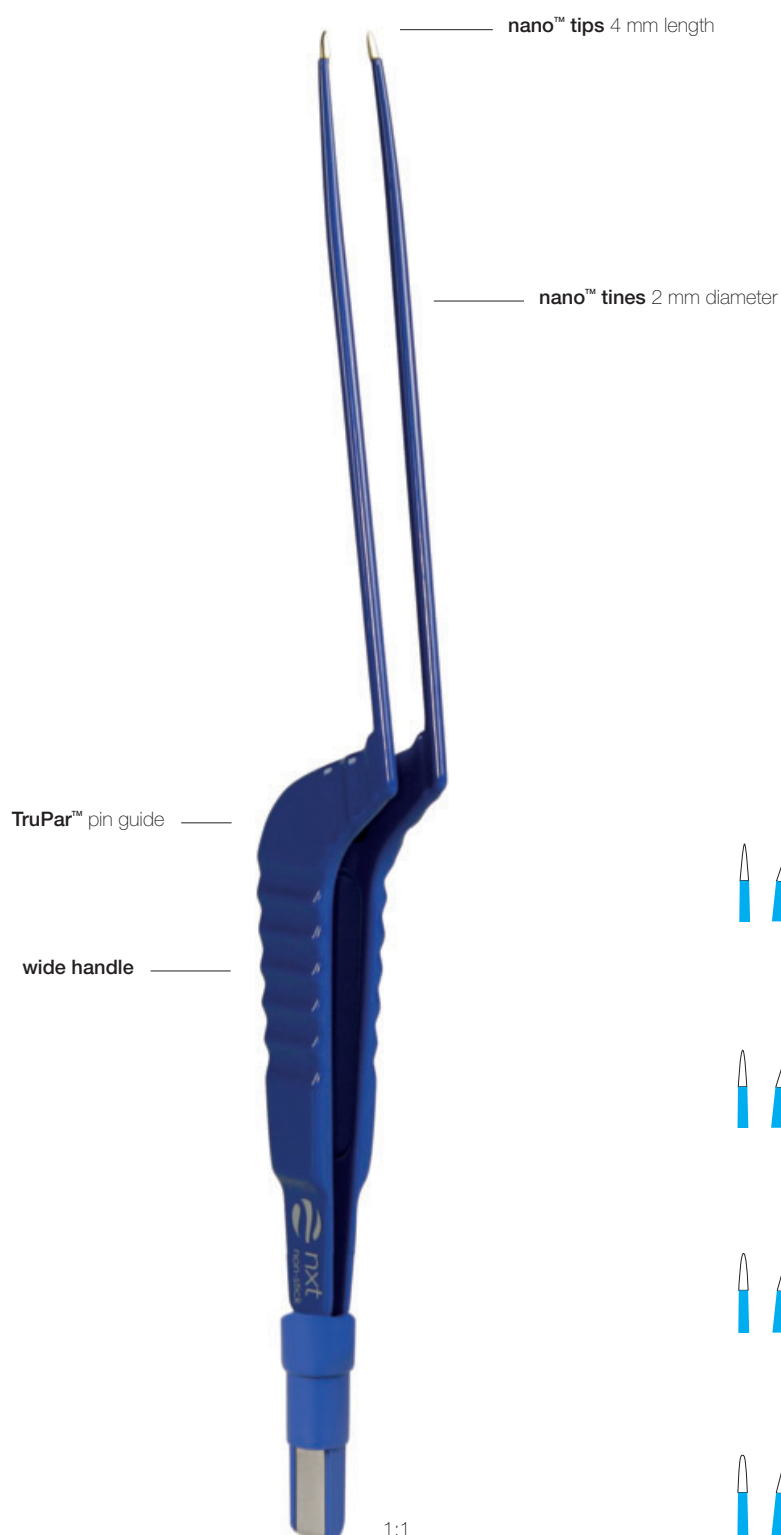

nano™ non-stick



Ø 0.2 x 4 mm nano™ tip straight

nano™ series / working / total length 85 / 200 mm

BIP004EU eu flat connector

BIP004US us pin connector



Ø 0.4 x 4 mm nano™ tip straight

nano™ series / working / total length 85 / 200 mm

BIP005EU eu flat connector

BIP005US us pin connector



Ø 0.7 x 4 mm nano™ tip straight

nano™ series / working / total length 85 / 200 mm

BIP006EU eu flat connector

BIP006US us pin connector




Ø 1.0 x 4 mm nano™ tip straight

nano™ series / working / total length 85 / 200 mm

BIP007EU eu flat connector

BIP007US us pin connector

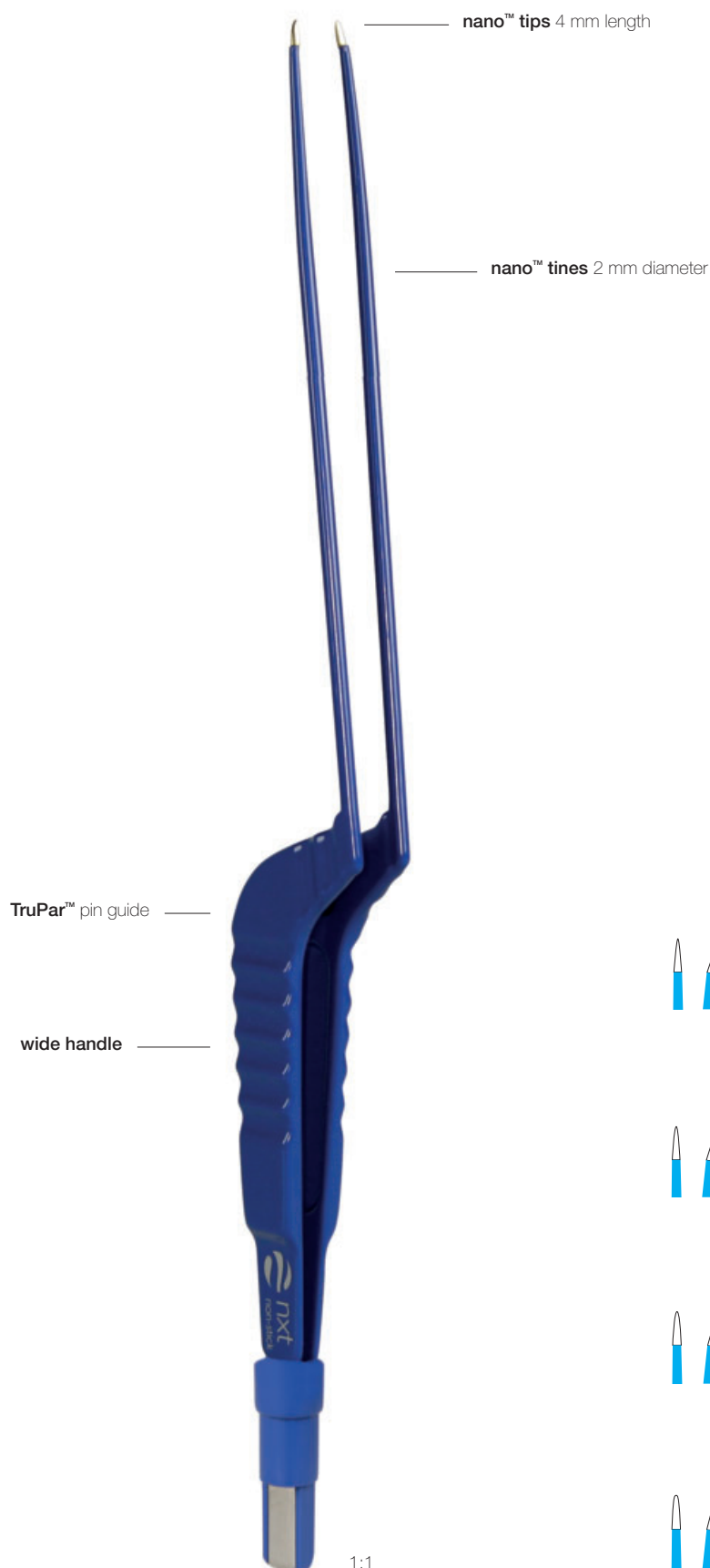

nano™ non-stick



Ø 0.2 x 4 mm nano™ tip straight

nano™ series / working / total length 115 / 230 mm

BIP008EU eu flat connector

BIP008US us pin connector



Ø 0.4 x 4 mm nano™ tip straight

nano™ series / working / total length 115 / 230 mm

BIP009EU eu flat connector

BIP009US us pin connector



Ø 0.7 x 4 mm nano™ tip straight

nano™ series / working / total length 115 / 230 mm

BIP010EU eu flat connector

BIP010US us pin connector



Ø 1.0 x 4 mm nano™ tip straight

nano™ series / working / total length 115 / 230 mm

BIP011EU eu flat connector

BIP011US us pin connector



nano™ non-stick



Ø 0.2 x 4 mm nano™ tip straight

nano™ series / working / total length 135 / 250 mm

BIP012EU eu flat connector

BIP012US us pin connector



Ø 0.4 x 4 mm nano™ tip straight

nano™ series / working / total length 135 / 250 mm

BIP013EU eu flat connector

BIP013US us pin connector



Ø 0.7 x 4 mm nano™ tip straight

nano™ series / working / total length 135 / 250 mm

BIP014EU eu flat connector

BIP014US us pin connector




Ø 1.0 x 4 mm nano™ tip straight

nano™ series / working / total length 135 / 250 mm

BIP015EU eu flat connector


BIP015US us pin connector

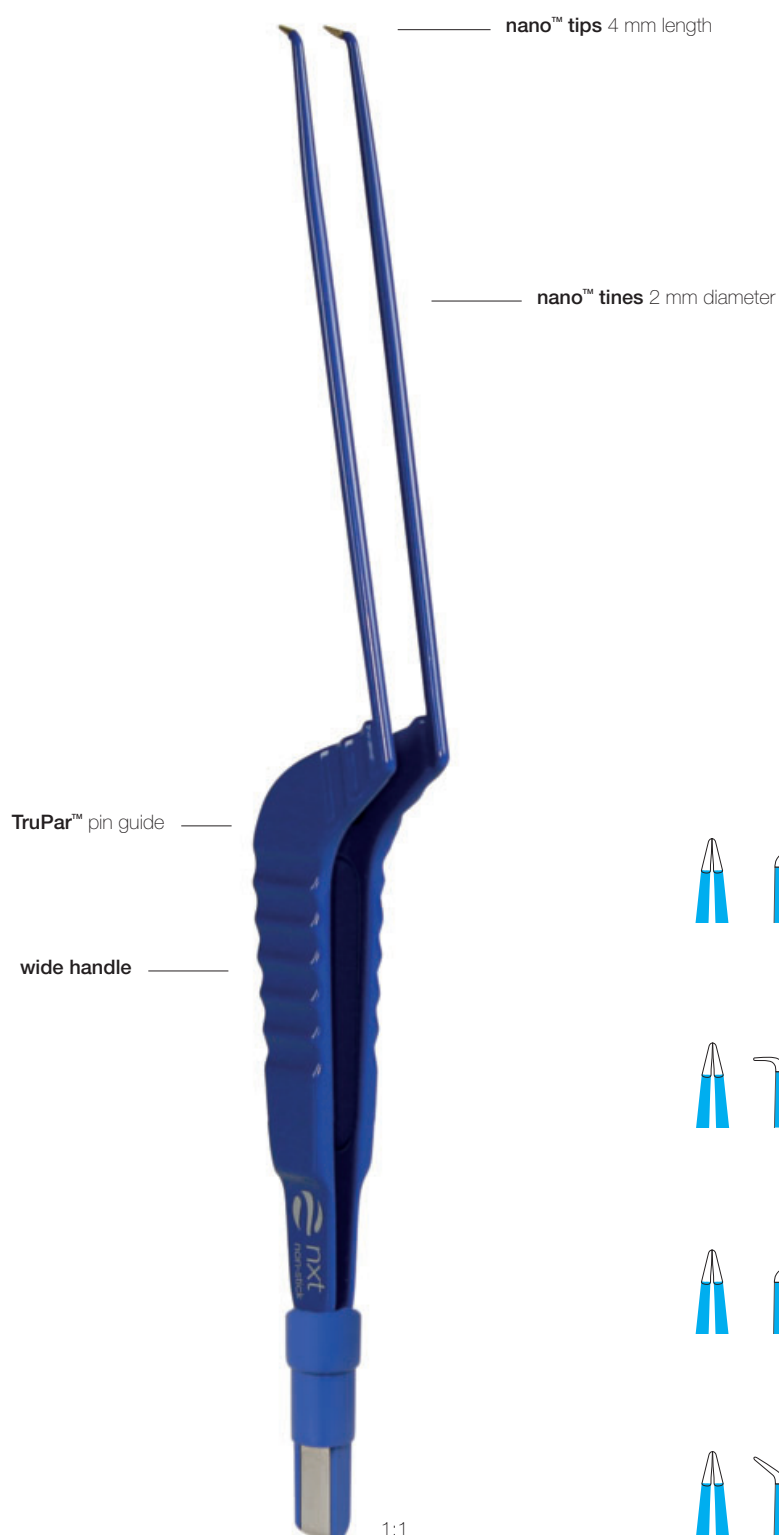

— nano™ tips 4 mm length

— nano™ tines 2 mm diameter


nano™ non-stick

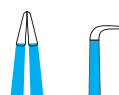

Ø 0.4 x 4 mm nano™ tip 80° upwards angled
nano™ series / working / total length 85 / 200 mm
BIP016EU eu flat connector
BIP016US us pin connector




Ø 0.4 x 4 mm nano™ tip 80° downwards angled
nano™ series / working / total length 85 / 200 mm
BIP017EU eu flat connector
BIP017US us pin connector

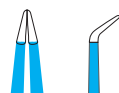


Ø 0.4 x 4 mm nano™ tip 45° upwards angled
nano™ series / working / total length 85 / 200 mm
BIP018EU eu flat connector
BIP018US us pin connector




Ø 0.4 x 4 mm nano™ tip 45° downwards angled
nano™ series / working / total length 85 / 200 mm
BIP019EU eu flat connector
BIP019US us pin connector

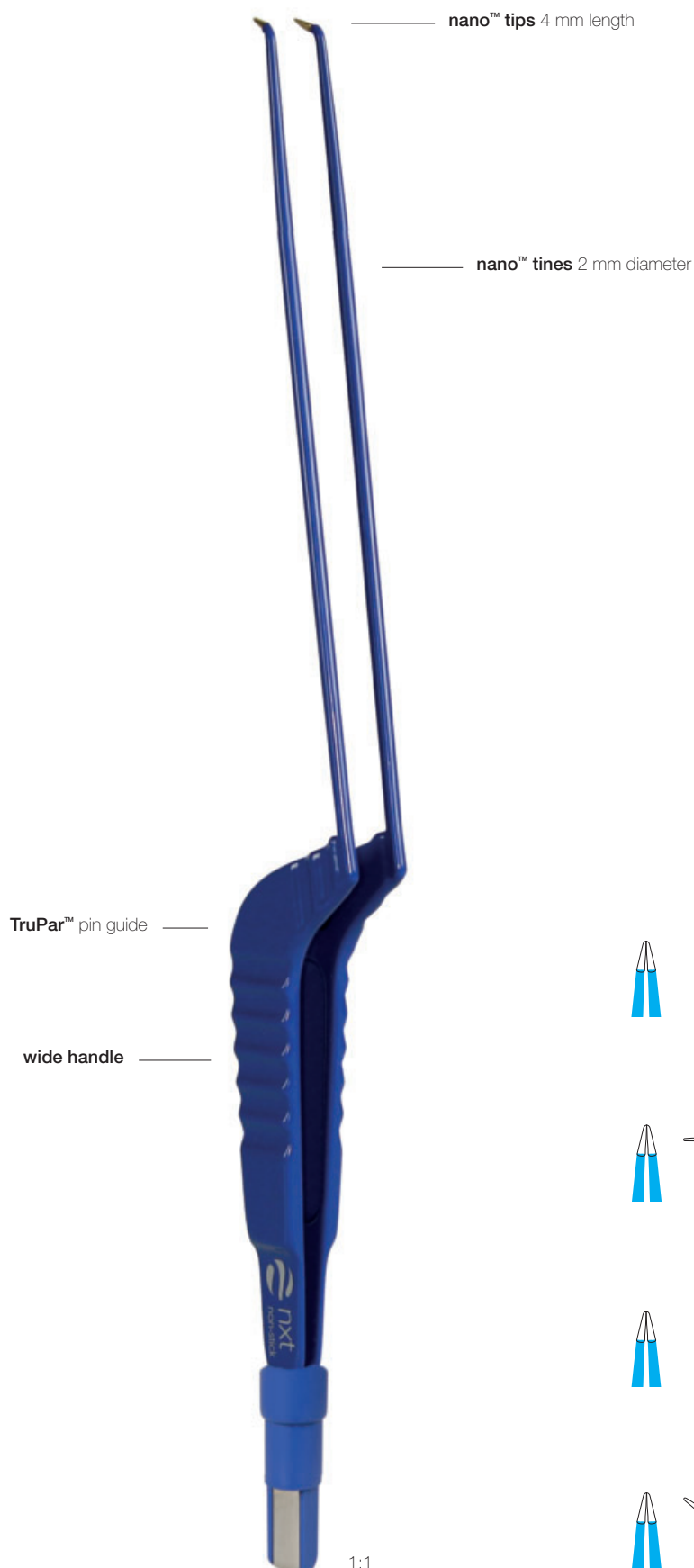


nano™ non-stick

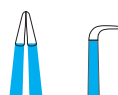


Ø 0.4 x 4 mm nano™ tip 80° upwards angled
nano™ series / working / total length 115 / 230 mm
BIP020EU eu flat connector
BIP020US us pin connector

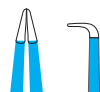

Ø 0.4 x 4 mm nano™ tip 80° downwards angled
nano™ series / working / total length 115 / 230 mm
BIP021EU eu flat connector
BIP021US us pin connector

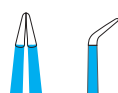

Ø 0.4 x 4 mm nano™ tip 45° upwards angled
nano™ series / working / total length 115 / 230 mm
BIP022EU eu flat connector
BIP022US us pin connector

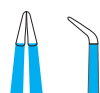


Ø 0.4 x 4 mm nano™ tip 45° downwards angled
nano™ series / working / total length 115 / 230 mm
BIP023EU eu flat connector
BIP023US us pin connector



1:1

classic non-stick

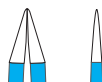


Ø 0.2 x 6 mm micro tip straight

standard series / working / total length 30 / 120 mm

BIP024EU eu flat connector

BIP024US us pin connector

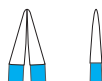


Ø 0.4 x 6 mm micro tip straight

standard series / working / total length 30 / 120 mm

BIP025EU eu flat connector

BIP025US us pin connector

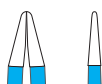


Ø 0.7 x 6 mm micro tip straight

standard series / working / total length 30 / 120 mm

BIP026EU eu flat connector

BIP026US us pin connector

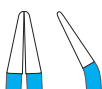


Ø 0.7 x 8 mm standard tip 45° angled

standard series / working / total length 30 / 120 mm

BIP027EU eu flat connector

BIP027US us pin connector

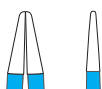

Ø 1.0 x 8 mm standard tip straight


standard series / working / total length 30 / 120 mm

BIP028EU eu flat connector

BIP028US us pin connector



classic non-stick

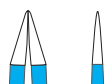


Ø 0.2 x 6 mm micro tip straight

standard series / working / total length 45 / 185 mm

BIP029EU eu flat connector

BIP029US us pin connector

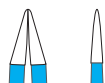


Ø 0.4 x 6 mm micro tip straight

standard series / working / total length 45 / 185 mm

BIP030EU eu flat connector

BIP030US us pin connector

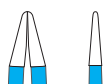


Ø 0.7 x 6 mm micro tip straight

standard series / working / total length 45 / 185 mm

BIP031EU eu flat connector

BIP031US us pin connector

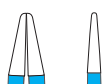

Ø 1.0 x 8 mm standard tip straight

standard series / working / total length 45 / 185 mm

BIP032EU eu flat connector

BIP032US us pin connector



Ø 1.0 x 8 mm standard tip 45° angled


standard series / working / total length 45 / 185 mm

BIP033EU eu flat connector

BIP033US us pin connector



classic non-stick

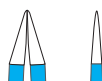


Ø 0.2 x 6 mm micro tip straight

standard series / working / total length 60 / 200 mm

BIP034EU eu flat connector

BIP034US us pin connector

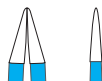


Ø 0.4 x 6 mm micro tip straight

standard series / working / total length 60 / 200 mm

BIP035EU eu flat connector

BIP035US us pin connector

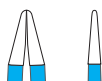


Ø 0.7 x 6 mm micro tip straight

standard series / working / total length 60 / 200 mm

BIP036EU eu flat connector

BIP036US us pin connector

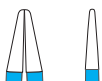

Ø 1.0 x 8 mm standard tip straight

standard series / working / total length 60 / 200 mm

BIP037EU eu flat connector

BIP037US us pin connector




Ø 1.0 x 8 mm standard tip 45° angled

standard series / working / total length 60 / 200 mm

BIP038EU eu flat connector

BIP038US us pin connector

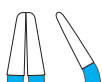

Ø 2.0 x 8 mm standard tip 45° angled

standard series / working / total length 60 / 200 mm

BIP039EU eu flat connector

BIP039US us pin connector

TruPar™ pin guide

wide handle



1:1

classic non-stick



Ø 1.5 x 8 mm standard tip straight

standard series / working / total length 90 / 230 mm

BIP040EU eu flat connector

BIP040US us pin connector



Ø 2.0 x 8 mm standard tip straight

standard series / working / total length 90 / 230 mm

BIP041EU eu flat connector

BIP041US us pin connector


TruPar™ pin guide

wide handle




1:1

classic non-stick

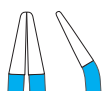

Ø 1.0 x 8 mm standard tip 45° angled
standard series / working / total length 120 / 260 mm
BIP042EU eu flat connector
BIP042US us pin connector

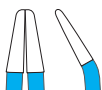

Ø 2.0 x 8 mm standard tip 45° angled
standard series / working / total length 120 / 260 mm
BIP043EU eu flat connector
BIP043US us pin connector

TruPar™ pin guide

wide handle



1:1

classic non-stick



Ø 2.5 x 8 mm standard tip straight


standard series / working / total length 140 / 280 mm

BIP044EU eu flat connector

BIP044US us pin connector



classic non-stick

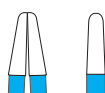

Ø 2.0 x 8 mm standard tip straight
standard series / working / total length 50 / 180 mm
BIP045EU eu flat connector
BIP045US us pin connector


TruPar™ pin guide

wide handle



1:1

classic non-stick

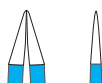


Ø 0.2 x 6 mm micro tip straight

standard series / working / total length 40 / 155 mm

BIP046EU eu flat connector

BIP046US us pin connector

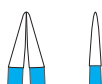


Ø 0.4 x 6 mm micro tip straight

standard series / working / total length 40 / 155 mm

BIP047EU eu flat connector

BIP047US us pin connector

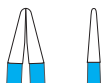


Ø 0.7 x 6 mm micro tip straight


standard series / working / total length 40 / 155 mm

BIP048EU eu flat connector


BIP048US us pin connector

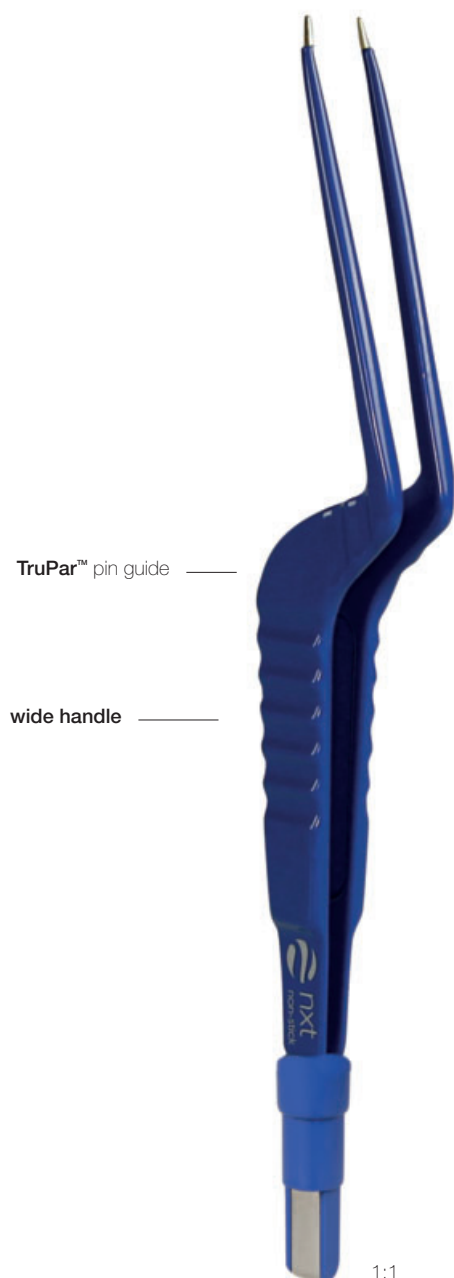


classic non-stick

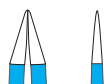


Ø 0.2 x 6 mm micro tip straight
standard series / working / total length 55 / 170 mm
BIP049EU eu flat connector
BIP049US us pin connector

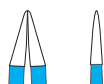


Ø 0.4 x 6 mm micro tip straight
standard series / working / total length 55 / 170 mm
BIP050EU eu flat connector
BIP050US us pin connector

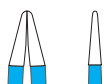

Ø 0.7 x 6 mm micro tip straight
standard series / working / total length 55 / 170 mm
BIP051EU eu flat connector
BIP051US us pin connector

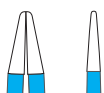

Ø 1.0 x 8 mm standard tip straight
standard series / working / total length 55 / 170 mm
BIP052EU eu flat connector
BIP052US us pin connector

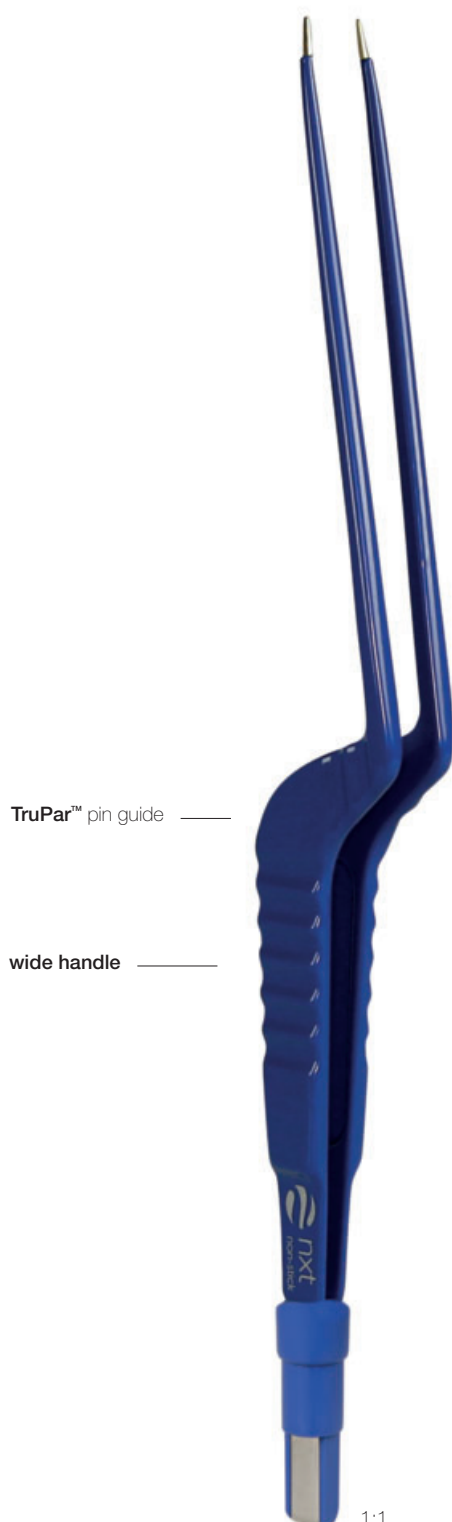

classic non-stick



Ø 0.2 x 6 mm micro tip straight

standard series / working / total length 85 / 200 mm

BIP053EU eu flat connector

BIP053US us pin connector



Ø 0.4 x 6 mm micro tip straight

standard series / working / total length 85 / 200 mm

BIP054EU eu flat connector

BIP054US us pin connector



Ø 0.7 x 6 mm micro tip straight

standard series / working / total length 85 / 200 mm

BIP055EU eu flat connector

BIP055US us pin connector



Ø 1.0 x 8 mm standard tip straight

standard series / working / total length 85 / 200 mm

BIP056EU eu flat connector

BIP056US us pin connector



Ø 1.2 x 8 mm standard tip straight

standard series / working / total length 85 / 200 mm

BIP057EU eu flat connector

BIP057US us pin connector




Ø 2.0 x 8 mm standard tip straight


standard series / working / total length 85 / 200 mm

BIP058EU eu flat connector


BIP058US us pin connector

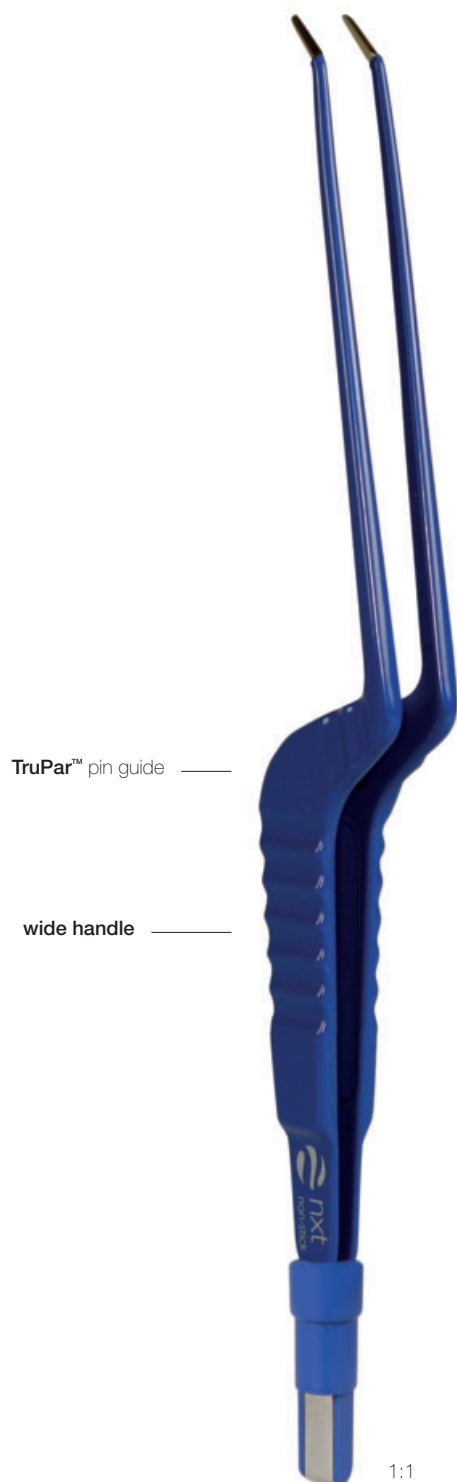

classic non-stick

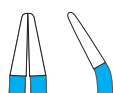

Ø 1.2 x 8 mm standard tip 45° downwards angled
standard series / working / total length 85 / 200 mm
BIP059EU eu flat connector
BIP059US us pin connector

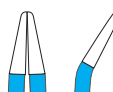

Ø 1.2 x 8 mm standard tip 45° upwards angled
standard series / working / total length 85 / 200 mm
BIP060EU eu flat connector
BIP060US us pin connector



classic non-stick



Ø 0.2 x 6 mm micro tip straight

standard series / working / total length 115 / 230 mm

BIP061EU eu flat connector

BIP061US us pin connector



Ø 0.4 x 6 mm micro tip straight

standard series / working / total length 115 / 230 mm

BIP062EU eu flat connector

BIP062US us pin connector



Ø 0.7 x 6 mm micro tip straight

standard series / working / total length 115 / 230 mm

BIP063EU eu flat connector

BIP063US us pin connector



Ø 1.0 x 8 mm standard tip straight

standard series / working / total length 115 / 230 mm

BIP064EU eu flat connector

BIP064US us pin connector



Ø 1.2 x 8 mm standard tip straight


standard series / working / total length 115 / 230 mm

BIP065EU eu flat connector


BIP065US us pin connector




classic non-stick

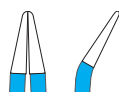


Ø 0.7 x 8 mm standard tip 45° angled upwards
standard series / working / total length 115 / 230 mm
BIP066EU eu flat connector
BIP066US us pin connector

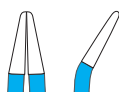

Ø 1.2 x 8 mm standard tip 45° angled upwards
standard series / working / total length 115 / 230 mm
BIP067EU eu flat connector
BIP067US us pin connector

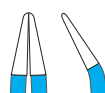

Ø 0.7 x 8 mm standard tip 45° angled downwards
standard series / working / total length 115 / 230 mm
BIP068EU eu flat connector
BIP068US us pin connector

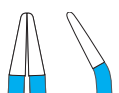

Ø 1.2 x 8 mm standard tip 45° angled downwards
standard series / working / total length 115 / 230 mm
BIP069EU eu flat connector
BIP069US us pin connector

TruPar™ pin guide

wide handle



1:1

classic non-stick



Ø 0.2 x 6 mm micro tip straight

standard series / working / total length 135 / 250 mm

BIP070EU eu flat connector

BIP070US us pin connector



Ø 0.4 x 6 mm micro tip straight

standard series / working / total length 135 / 250 mm

BIP071EU eu flat connector

BIP071US us pin connector



Ø 0.7 x 6 mm micro tip straight

standard series / working / total length 135 / 250 mm

BIP072EU eu flat connector

BIP072US us pin connector



Ø 1.0 x 8 mm standard tip straight

standard series / working / total length 135 / 250 mm

BIP073EU eu flat connector

BIP073US us pin connector



Ø 1.2 x 8 mm standard tip straight

standard series / working / total length 135 / 250 mm

BIP074EU eu flat connector

BIP074US us pin connector


TruPar™ pin guide

wide handle




1:1


classic non-stick

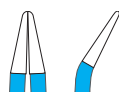


Ø 0.7 x 8 mm standard tip 45° angled upwards
standard series / working / total length 135 / 250 mm
BIP075EU eu flat connector
BIP075US us pin connector

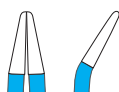

Ø 1.2 x 8 mm standard tip 45° angled upwards
standard series / working / total length 135 / 250 mm
BIP076EU eu flat connector
BIP076US us pin connector

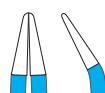

Ø 0.7 x 8 mm standard tip 45° angled downwards
standard series / working / total length 135 / 250 mm
BIP077EU eu flat connector
BIP077US us pin connector

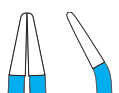

Ø 1.2 x 8 mm standard tip 45° angled downwards
standard series / working / total length 135 / 250 mm
BIP078EU eu flat connector
BIP078US us pin connector



non-stick with integrated irrigation




Integrated irrigation channel
 Integrierter Spülkanal




TruPar™ Pin guide




Standardized irrigation connector
 Standardisierter Spülanschluss

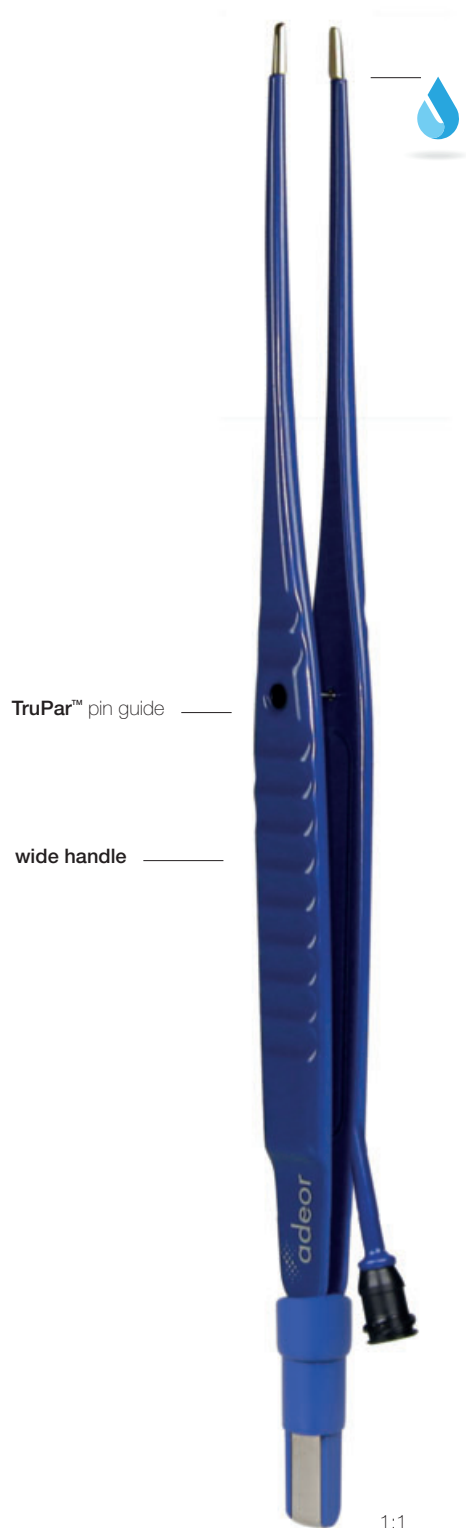


irrigation non-stick

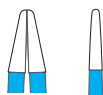


Ø 1.5 x 8 mm standard tip straight
standard series / working / total length 100 / 200 mm
BIP001EU-IR eu flat connector
BIP001US-IR us pin connector

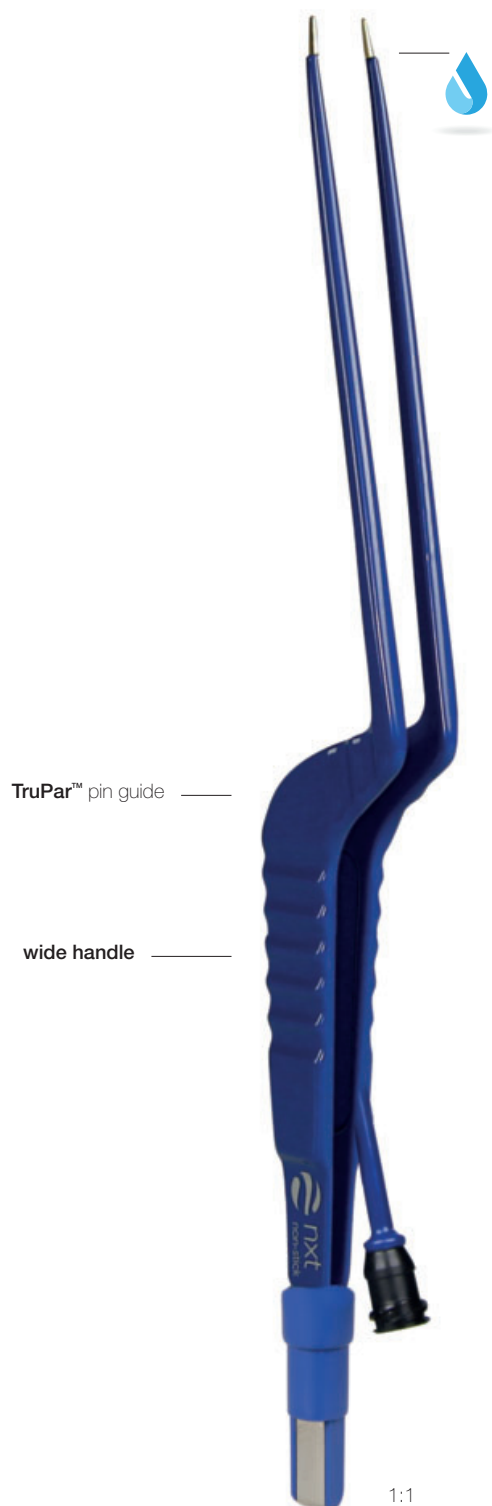

irrigation non-stick

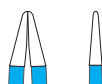


Ø 0.5 x 6 mm micro tip straight
standard series / working / total length 90 / 200 mm
BIP002EU-IR eu flat connector
BIP002US-IR us pin connector

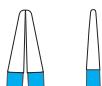


Ø 1.0 x 8 mm standard tip straight
standard series / working / total length 90 / 200 mm
BIP003EU-IR eu flat connector
BIP003US-IR us pin connector



irrigation non-stick

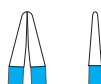

Ø 0.5 x 6 mm micro tip straight
standard series / working / total length 120 / 230 mm
BIP004EU-IR eu flat connector
BIP004US-IR us pin connector

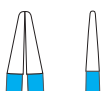

Ø 1.0 x 8 mm standard tip straight
standard series / working / total length 120 / 230 mm
BIP005EU-IR eu flat connector
BIP005US-IR us pin connector



suction forceps



Integrated suction channel
 Integrierter Saugkanal



Wide ergonomic grip
 Breiter Ergo-Griff



Standardized suction connector
 Standardisierter Sauganschluss


TruPar™ pin guide

wide handle




1:1

suction forceps

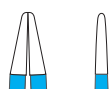


Ø 1.4 mm standard tip straight
standard series / working / total length 90 / 200 mm
BIP001EU-SUC eu flat connector
BIP001US-SUC us pin connector

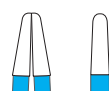

Ø 2.0 mm standard tip straight
standard series / working / total length 90 / 200 mm
BIP002EU-SUC eu flat connector
BIP002US-SUC us pin connector



suction forceps

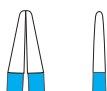


Ø 1.4 mm standard tip straight

standard series / working / total length 120 / 230 mm

BIP003EU-SUC eu flat connector

BIP003US-SUC us pin connector

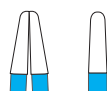

Ø 2.0 mm standard tip straight


standard series / working / total length 120 / 230 mm

BIP004EU-SUC eu flat connector


BIP004US-SUC us pin connector



suction forceps

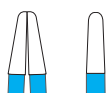

Ø 2.0 standard tip straight
standard series / working / total length 75 / 190 mm
BIP005EU-SUC eu flat connector
BIP005US-SUC us pin connector

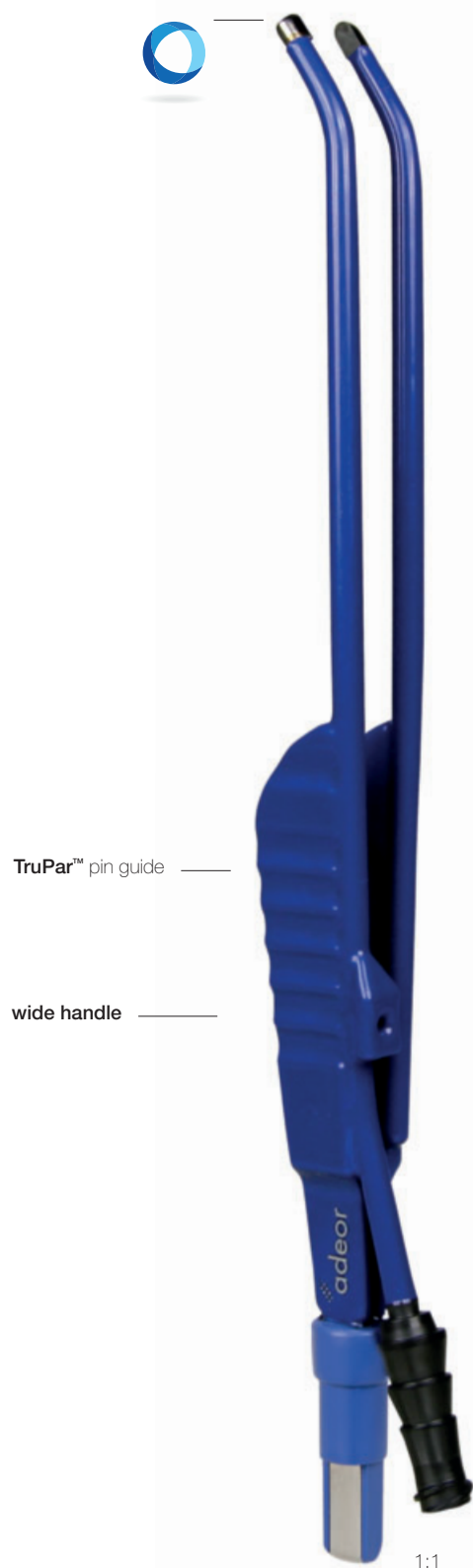

meuser suction forceps



Ø 3.0 standard tip straight


standard series / working / total length 75 / 210 mm

BIP006EU-SUC eu flat connector

BIP006US-SUC us pin connector



meuser suction forceps

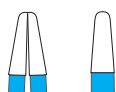


Ø 3.2 standard tip straight
standard series / working / total length 75 / 210 mm
BIP007EU-SUC eu flat connector
BIP007US-SUC us pin connector




ZEPELIN® style round-handle non-stick



nano™ tips (4 mm length)
 nano™ Spitzen (4 mm Länge)

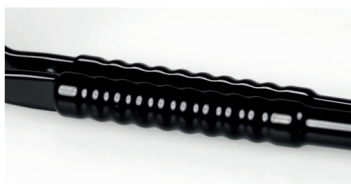

Round-handle design
 Rundgriff-Design

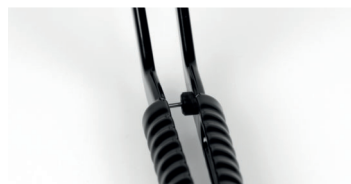

TruPar™ Pin guide



Completely closed Hygiene+™ connector
 Geschlossener Hygiene+™ Stecker



Ø 0.5 x 4 mm straight

working / total length 50 / 150 mm

BIP-S141EU eu flat connector

BIP-S141US us pin connector



Ø 0.8 x 4 mm straight

working / total length 50 / 150 mm

BIP-S142EU eu flat connector

BIP-S142US us pin connector



Ø 0.5 x 4 mm angled

working / total length 50 / 150 mm

BIP-S146EU eu flat connector

BIP-S146US us pin connector



Ø 0.8 x 4 mm angled

working / total length 50 / 150 mm

BIP-S147EU eu flat connector

BIP-S147US us pin connector



Ø 0.5 x 4 mm angled

working / total length 80 / 180 mm

BIP-S246EU eu flat connector

BIP-S246US us pin connector



Ø 1.2 x 4 mm angled

working / total length 80 / 180 mm

BIP-S248EU eu flat connector

BIP-S248US us pin connector



Ø 0.5 x 4 mm straight

working / total length 100 / 200 mm

BIP-S341EU eu flat connector

BIP-S341US us pin connector



Ø 0.8 x 4 mm straight

working / total length 100 / 200 mm

BIP-S342EU eu flat connector

BIP-S342US us pin connector



Ø 0.5 x 4 mm angled

working / total length 100 / 200 mm

BIP-S346EU eu flat connector

BIP-S346US us pin connector




Ø 0.8 x 4 mm angled

working / total length 100 / 200 mm

BIP-S347EU eu flat connector

BIP-S347US us pin connector

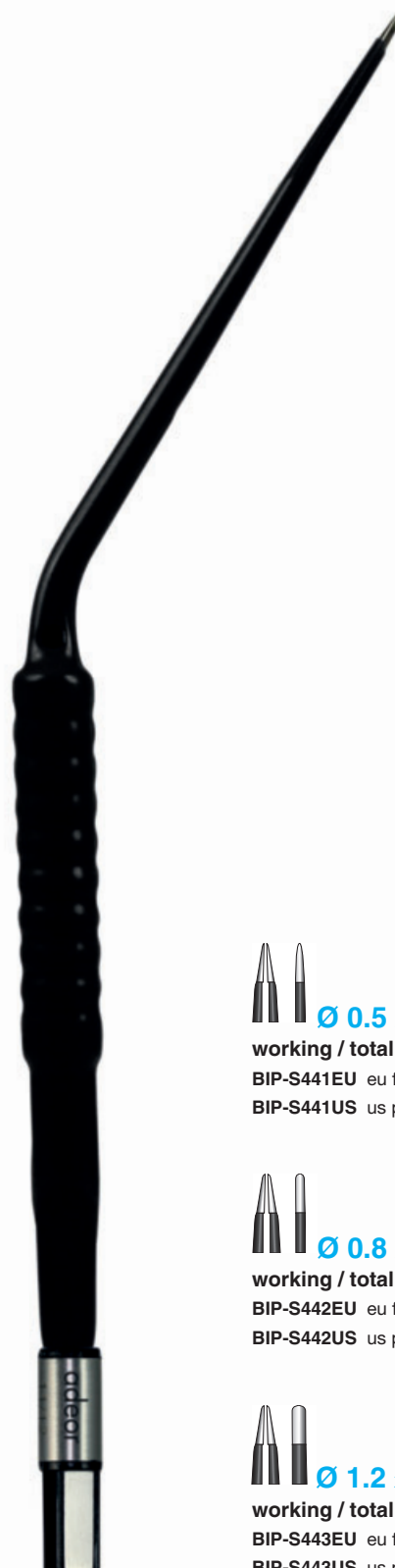

Ø 0.5 x 4 mm straight

working / total length 90 / 200 mm

BIP-S441EU eu flat connector

BIP-S441US us pin connector



Ø 0.8 x 4 mm straight

working / total length 90 / 200 mm

BIP-S442EU eu flat connector

BIP-S442US us pin connector



Ø 1.2 x 4 mm straight

working / total length 90 / 200 mm

BIP-S443EU eu flat connector

BIP-S443US us pin connector



Ø 0.5 x 4 mm down angled
working / total length 100 / 195 mm
BIP-S444EU eu flat connector
BIP-S444US us pin connector



Ø 0.5 x 4 mm up angled
working / total length 90 / 200 mm
BIP-S446EU eu flat connector
BIP-S446US us pin connector



Ø 0.8 x 4 mm up angled
working / total length 90 / 200 mm
BIP-S447EU eu flat connector
BIP-S447US us pin connector



Ø 1.2 x 4 mm up angled
working / total length 90 / 200 mm
BIP-S448EU eu flat connector
BIP-S448US us pin connector



Ø 0.5 x 4 mm straight

working / total length 140 / 240 mm

BIP-S541EU eu flat connector

BIP-S541US us pin connector



Ø 0.8 x 4 mm straight

working / total length 140 / 240 mm

BIP-S542EU eu flat connector

BIP-S542US us pin connector



Ø 0.5 x 4 mm angled

working / total length 140 / 240 mm

BIP-S546EU eu flat connector

BIP-S546US us pin connector



Ø 0.8 x 4 mm angled

working / total length 140 / 240 mm

BIP-S547EU eu flat connector

BIP-S547US us pin connector



Ø 0.2 x 4 mm straight

working / total length 75 / 220 mm

BIP-S641EU eu flat connector

BIP-S641US us pin connector



Ø 0.5 x 4 mm straight

working / total length 75 / 220 mm

BIP-S642EU eu flat connector

BIP-S642US us pin connector



Ø 0.8 x 4 mm straight

working / total length 75 / 220 mm

BIP-S643EU eu flat connector

BIP-S643US us pin connector



Ø 1.2 x 4 mm straight

working / total length 75 / 220 mm

BIP-S644EU eu flat connector

BIP-S644US us pin connector



Ø 2.0 x 4 mm straight

working / total length 75 / 220 mm

BIP-S645EU eu flat connector

BIP-S645US us pin connector



Ø 0.5 x 4 mm up angled

working / total length 80 / 220 mm

BIP-S646EU eu flat connector

BIP-S646US us pin connector



Ø 0.8 x 4 mm up angled

working / total length 80 / 220 mm

BIP-S647EU eu flat connector

BIP-S647US us pin connector



Ø 1.2 x 4 mm up angled

working / total length 80 / 220 mm

BIP-S648EU eu flat connector

BIP-S648US us pin connector



Ø 0.5 x 4 mm straight

working / total length 110 / 250 mm

BIP-S741EU eu flat connector

BIP-S741US us pin connector



Ø 0.8 x 4 mm straight

working / total length 110 / 250 mm

BIP-S742EU eu flat connector

BIP-S742US us pin connector



Ø 1.2 x 4 mm straight

working / total length 110 / 250 mm

BIP-S743EU eu flat connector

BIP-S743US us pin connector



Ø 0.5 x 4 mm up angled

working / total length 110 / 250 mm

BIP-S746EU eu flat connector

BIP-S746US us pin connector



Ø 0.8 x 4 mm up angled

working / total length 110 / 250 mm

BIP-S747EU eu flat connector

BIP-S747US us pin connector



Ø 1.2 x 4 mm up angled


working / total length 110 / 250 mm

BIP-S748EU eu flat connector

BIP-S748US us pin connector

BIPOLAR CABLES

RE-USABLE NON-STERILE

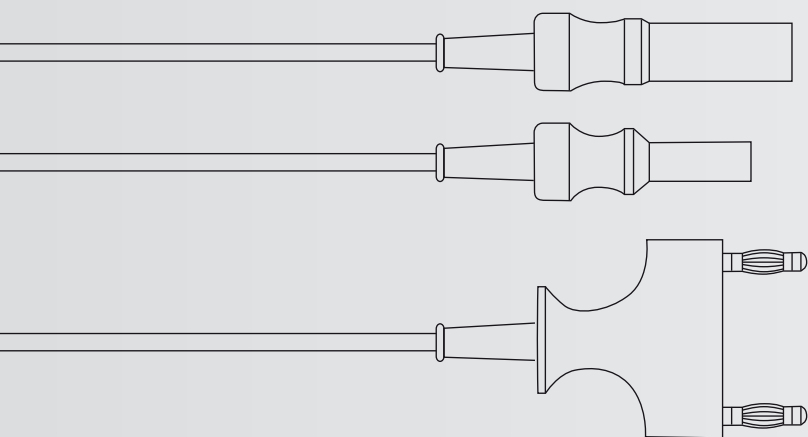

bipolar cables*

*not an adeor product.

*kein adeor Produkt.

Generator side



Forceps side

	EU flat	US pin
Valleylab, Conmed, Codman, Olympus - Banana plug*	ST-92021 3 m	ST-92114 3 m
Valleylab, Conmed, Codman, Olympus - Banana plug*	ST-92022 5 m	ST-92093 5 m
Aesculap, Martin, Berchtold, Wolf	ST-92017 3 m	ST-92113 3 m
Aesculap, Martin, Berchtold, Wolf	ST-92018 5 m	ST-92092 5 m
Erbe, Downs, Siemens	ST-92015 3 m	ST-92110 3 m
Erbe, Downs, Siemens	ST-92016 5 m	ST-92111 5 m
Valleylab, EMD, Doley - 2-pin distance 29 mm	ST-92019 3 m	ST-92115 3 m
Valleylab, EMD, Doley - 2-pin distance 29 mm	ST-92020 5 m	ST-92094 5 m

② multi use
autoclavable 134°C

*not available in the EU and some other markets.

*nicht in der EU und einigen anderen Märkten verfügbar.

HEADQUARTERS

adeor medical AG
Martinshof 5
D-83626 Valley
Germany

Phone +49 (0) 8024 47717 0
Fax +49 (0) 8024 47717 99
eMail office@adeor.com
web www.adeor.com



Printed on recycled paper.

Gedruckt auf Recyclingpapier.

